

GCC Lab المختبر الفليبي	Gulf Laboratory Technical Services
	Operational Procedure of Color meter Lovi Bond EC 3000ASTM

Equipment Number: EQ-05

Serial Number: 201360



1. PURPOSE

This standard operating procedure is intended to provide users with the guide on the operation and maintenance of **Color meter Lovi Bond EC 3000ASTM**.

2. RESPONSIBILITIES

- Use, care and maintenance of the **Color meter Lovi Bond EC 3000ASTM** is directly assigned to the end user who is usually the laboratory chemist. This responsibility may be shared to other qualified user(s) upon the consent of the laboratory chemist.

3. PROCEDURE

EC 3000 ASTM Color (ASTM D1500): because in the scope of ASTM: D6045 it states, “This test method correlates to Test Method D 156 and Test Method D 1500 as calculated by the instrumentation”.

Parameter	Details
Colour Scales	ASTM
Range	0 to 8
Resolution	0.1
Repeatability	± 0.1

Doc.: EOP [6]-04-F	Version No.00	Issue Date: 1/11/2023	Revision No: 00	Date:1/11/2023
Prepared by: Quality Manager	Reviewed by: Technical Manager	Approved by: General Manager	Page 1 of 4	

GCC Lab المفتبر الفليبي	Gulf Laboratory Technical Services
	Operational Procedure of Color meter Lovi Bond EC 3000ASTM

Performing a Zero - ASTM:




- After initiation, the instrument will display:
- Ensure the sample chamber is empty.
- Close the sample lid and press the zero key on the keypad or the display.
- While performing a zero, the instrument will display.
- When complete, the instrument will display.



Performing a Test - ASTM:

Note: Ensure that the sample being measured is free from turbidity; if necessary, filter the sample through a 0.45 micron filter paper.


- Fill a clean Cell 2/3 full with the sample.
- Place the cell in the sample chamber of the EComparator, push to the left hand side.
- Close the sample lid and press the test key on the keypad or the display.
- While performing a test, the instrument will display.
- The instrument will display the result as
- Press the result button on the display to change between numeric and comparator mode.

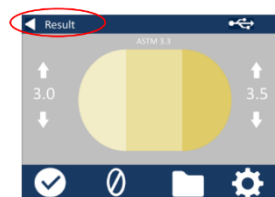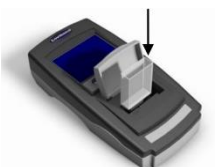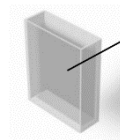

Doc.: EOP [6]-04-F	Version No.00	Issue Date: 1/11/2023	Revision No: 00	Date:1/11/2023
Prepared by: Quality Manager	Reviewed by: Technical Manager	Approved by: General Manager	Page 2 of 4	

GCC Lab المفتبر الفليبي	Gulf Laboratory Technical Services
	Operational Procedure of Color meter Lovi Bond EC 3000ASTM

- Press the up and down arrows to view the stages / steps as they would be on the visual disc.

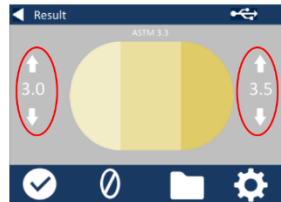


Performing a Liquid Instrument Validation using a Liquid Reference Standard – ASTM:

The EC 3000 ASTM is supplied with two bottles of Validation Liquid.

These liquids are used to validate that the instrument is measuring within the specified tolerances.



- Select a bottle of validation liquid and ensure that the selected validation liquid is within its expiry date.

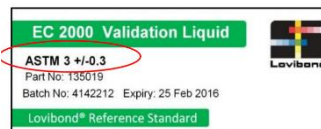


- Ensure the sample chamber is empty.
- Close the sample lid and press the zero key on the keypad or the display.
- While performing a zero, the instrument will display.



- When complete, the instrument will display.



- Fill a clean Cell 2/3 full with the reference solution.
- Place the cell in the sample chamber of the EComparator, push to the left hand side.
- Close the sample lid and press the test key on the keypad or the display.

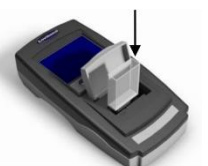

- While performing a test, the instrument will display.




Doc.: EOP [6]-04-F	Version No.00	Issue Date: 1/11/2023	Revision No: 00	Date:1/11/2023
Prepared by: Quality Manager	Reviewed by: Technical Manager	Approved by: General Manager	Page 3 of 4	

GCC Lab المفتبر الفليبي	Gulf Laboratory Technical Services
	Operational Procedure of Color meter Lovi Bond EC 3000ASTM

- The instrument will display the result as




- Compare the result from the instrument to the value on the label of the validation liquid, to check it is within the specified tolerances.

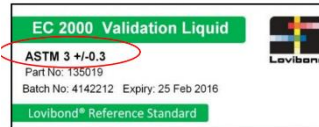

Performing an Instrument Validation using a Glass Reference Standard – ASTM

- Ensure the sample chamber is empty.

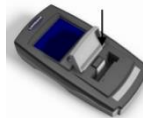

- Close the sample lid and press the zero key on the keypad or the display.



- Place the selected standard in the sample chamber on the right-hand side.



- Press the test key on the keypad.



- The result will be displayed on the screen. Compare these results to those on the label of the standard.



- Press the back arrow on the display to change between numeric and comparator mode



4. REFERENCES

[1] Operator Instruction Manual

Doc.: EOP [6]-04-F	Version No.00	Issue Date: 1/11/2023	Revision No: 00	Date:1/11/2023
Prepared by: Quality Manager	Reviewed by: Technical Manager	Approved by: General Manager	Page 4 of 4	